

Kul-PCR protokoll

Datum: 081003

Projekt: Alcedo / Neutrodenmatida

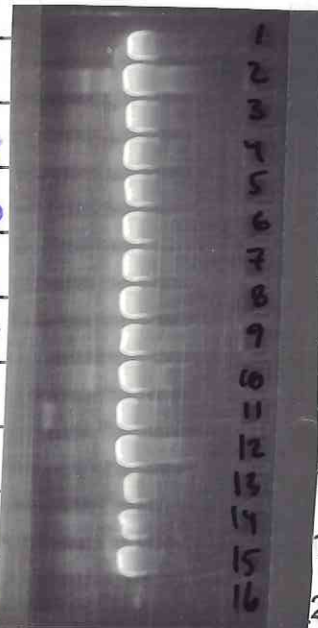
Primers: 1200F / 2450R

	µl/prov	mix x
Vatten	22	
Primer 1	0.15	
Primer 2	0.15	
BSA		
DNA	2	

Program:	File
95°C 3'	
94°C 30"	
58°C 30"	X40
72°C 1'	
72°C 10'	

Maskin nr: 23110 P

Taxa:

1. <u>08107</u>		11. <u>08119</u>
2. <u>08108</u>		12. <u>08120</u>
3. <u>08109</u>		13. <u>08121</u>
4. <u>08110</u>		14. <u>08122</u>
5. <u>08111</u>		15. <u>08123</u>
6. <u>08113</u>		16. <u>08126</u>
7. <u>08115</u>		17. <u>08131</u>
8. <u>08116</u>		18. <u>08-119 alt</u>
9. <u>08117</u>		19. <u>Kontroll</u>
10. <u>08118</u>		20. _____

Kommentarer: