

PCR protokoll

Datum: 111114 Projekt: Acoela-hösten2011

Primers: S30/5FK-acela - 4FB-acela/1806R


	ul/prov	mix X
H ₂ O	22	
primer 1	0,5	
primer 2	0,5	
DNA	2	

Program:

95°C	4min	X1
94°C	30"	X2
58-54°C	30"	
72°C	30"	
94°C	30"	x34
52°C	30"	
72°C	30"	
72°C	10 min	X1

Maskin nr:

1	UJ10-212
2	UJ10-221
3	UJ10-249
4	UJ11-047
5	UJ11-089
6	UJ11-318
7	UJ11-379
8	Kontroll
9	
10	
11	
12	
13	
14	
15	
16	
17	
18	
19	
20	
21	
22	
23	
24	

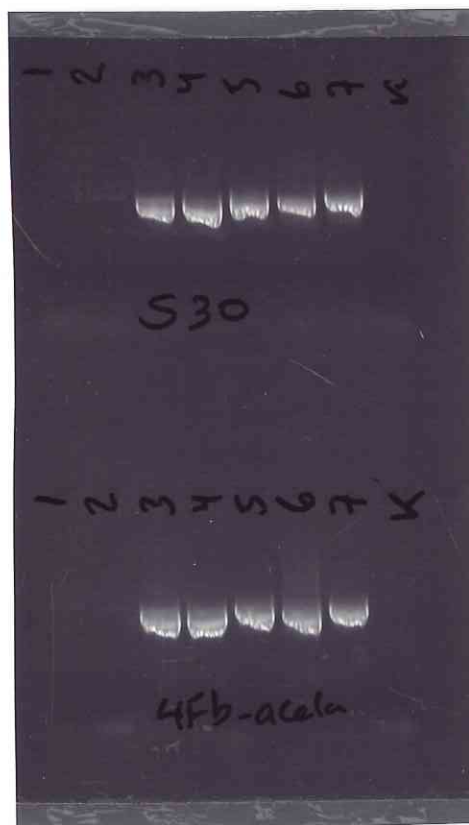

Ta hänsyn till att;

Dessa nya primers är något art specifika.

/Kerwan.