

PCR protokoll

Datum: 120217

Projekt: Acoela 2011

Primers: UJ2176//3449L


Program:

	ul/prov	mix X
H ₂ O	22	
primer 1	0,5	
primer 2	0,5	
DNA	2	

95°C 4min	X1
94°C 30"	X2
58-54 °C 30"	
72°C 1 min	
94C 30"	X34
52C 30"	
72C 1 min	
72°C 10 min	X1

Maskin nr:

1	UJ11_01	25	222	gelren
2	UJ11_17	26	224	gelren
3	UJ11_38	27	228	
4	UJ11_47	28	231	
5	UJ11_67	29	232	gelren
6	UJ11_74	30	234	gelren
7	UJ11_84	31	237	
8	UJ11_87	32	Kontroll	
9	UJ11_89			
10	UJ11_98			
11	UJ11_118			
12	UJ11_134			
13	UJ11_137			
14	UJ11_145			gelren
15	UJ11_149			
16	UJ11_169			
17	UJ11_179			
18	UJ11_191			
19	UJ11_201			
20	UJ11_203			
21	UJ11_204			
22	UJ11_205			
23	UJ11_220			
24	UJ11_221			

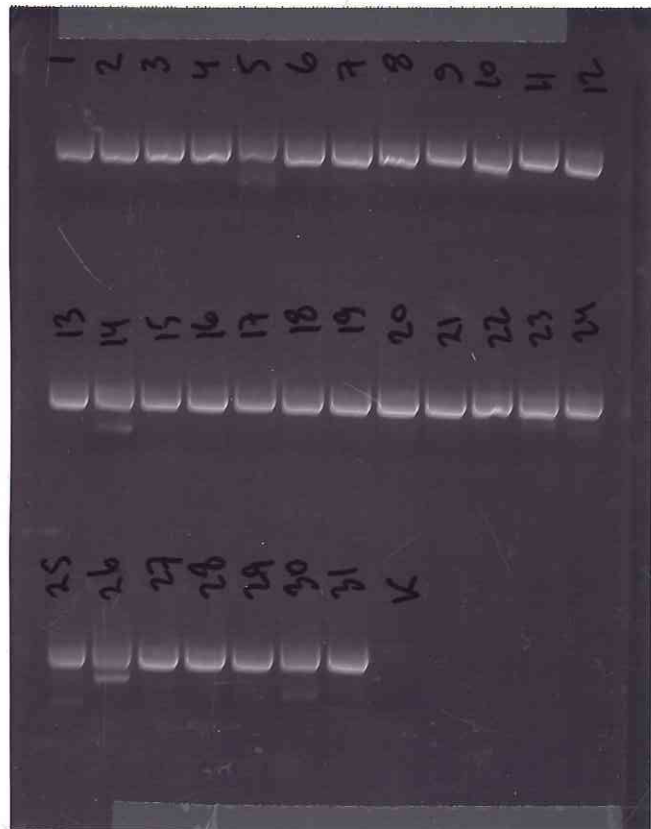

Keyvan