

PCR protokoll

Datum: 120222

Projekt: Acoela 2011

Primers: UJ2176//3449L


	ul/prov	mix X
H ₂ O	22	
primer 1	0,5	
primer 2	0,5	
DNA	2	

Program:

95°C 4min	X1
94°C 30"	X2
58-54 °C 30"	
72°C 1 min	
94C 30"	X34
52C 30"	
72C 1 min	
72°C 10 min	X1

Maskin nr:

1	UJ11_246	25	UJ11_344
2	UJ11_251	26	UJ11_356 <i>gelreana</i>
3	UJ11_252	27	UJ11_357
4	UJ11_258	28	UJ11_362
5	UJ11_259	29	UJ11_379
6	UJ11_260	30	UJ11_382 <i>gelreana</i>
7	UJ11_263 <i>gelreana</i>	31	UJ11_388
8	UJ11_272	32	UJ11_393
9	UJ11_302 <i>gelreana</i>	33	Kontroll
10	UJ11_303 <i>gelreana</i>		
11	UJ11_305		
12	UJ11_309 <i>gelreana</i>		
13	UJ11_314		
14	UJ11_315		
15	UJ11_319		
16	UJ11_320		
17	UJ11_324		
18	UJ11_325		
19	UJ11_329		
20	UJ11_330 <i>gelreana</i>		
21	UJ11_331		
22	UJ11_332		
23	UJ11_333		
24	UJ11_341 <i>gelreana</i>		

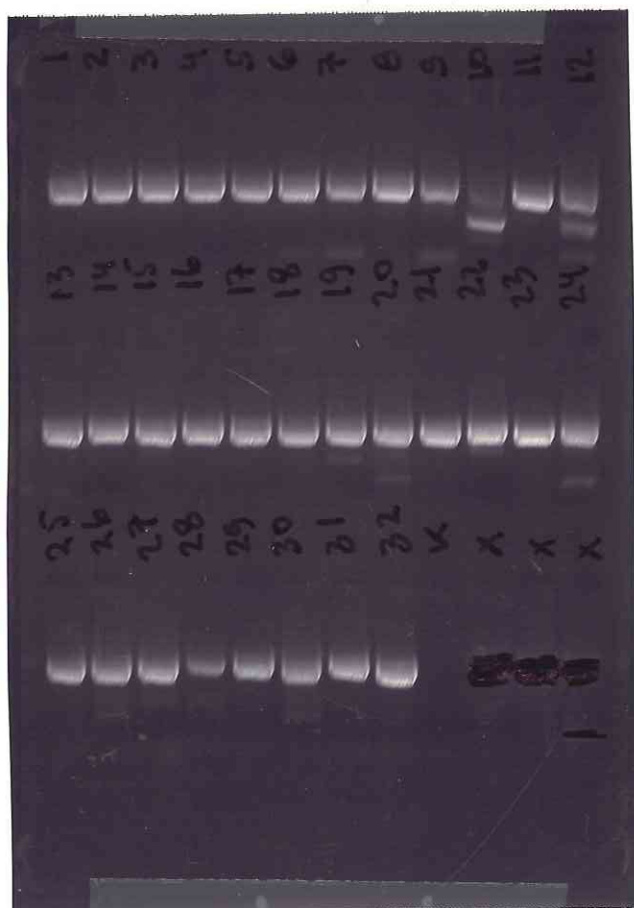

Kontroll:

Keyvan